



# Dental Nitrous Oxide Sedation




**MADE  
IN  
GERMANY**



# WELCOME

to Baldus Medizintechnik GmbH

» Your certified Specialist for Nitrous  
Oxide Sedation and Medical Gases

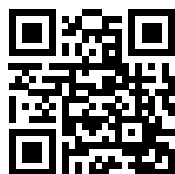

**Baldus Medizintechnik GmbH**

Auf dem Schafstall 5 · 56182 Urbar – Germany

☎ +49 (0) 261 / 96 38 926 - 0

☎ +49 (0) 261 / 96 38 926 - 21

✉ info@baldus-sedation.de

[www.baldus-sedation.de](http://www.baldus-sedation.de)

# CONTENTS

	<b>Dental Nitrous Oxide Sedation</b>	<b>4</b>
	Oxygen-Nitrous Oxide Mixer	6
	Comprehensive Set for Nitrous Oxide Sedation	12
	Baldus® Scavenger Breathing System	14
	Baldus® Single Use Inner Masks	20
	Baldus® Mounting Systems	23
	Sedation Equipment	26
	Scavenging Circuit Components	32
	<b>Emergency Equipment</b>	<b>35</b>
	<b>Central Gas Supply</b>	<b>37</b>
	<b>Maintenance</b>	<b>40</b>
	<b>n<sub>2</sub>o-Academy Workshops</b>	<b>42</b>

O<sub>2</sub> Flush

35 % N<sub>2</sub>O

7.0 l/min Gesamt-Flow

Bedienung auslösen

Menu

CE 0483 Made in Germany

Do. 12. Jan. 2017 13:10:10 Uhr

Behandlungszeit: 00h:10min 03sek

Medizintechnik Baldus Touch





# Dental Nitrous Oxide Sedation

## Anxiety free experience with Nitrous Oxide Sedation

The German Society of Oral and Maxillofacial Medicine estimates that up to 61% of the population fear the visit to the dentist. 20% of this percentage are considered to be highly anxious, and up to 10% avoid dentist completely due to severe anxiety.\*

Baldus has the solution for you and your patients:

Inhalation sedation with nitrous oxide and oxygen – is suitable for all ages. Your patient's fear or anxiety will be greatly reduced which will make for a far better atmosphere for both the patient and dentist. This will also allow for a better and more efficient treatment.

Nitrous oxide, or  $N_2O$ , has been used successfully for more than 160 years in dentistry. The mix of  $N_2O$  and  $O_2$  is the safest most researched and oldest sedation method. It's safe and non-toxic and when treatment is finished, the patient is given 100% oxygen, and within 15 minutes after session the patient is  $N_2O$  free. Compared to general anesthesia, the patient is fully conscious and able to communicate.


How does nitrous oxide work? During the treatment oxygen and nitrous oxide are mixed which is individually to the patient's needs. It starts with 100% oxygen then  $N_2O$  is raised successively until desired sedation is reached. While the patient is inhaling and exhaling the oxygen-nitrous oxide mixture, you and your team can start the treatment. Finally, the patient inhales 100% oxygen again, and within 15 minutes after session the patient is  $N_2O$  free.

## The advantages of Nitrous Oxide Sedation

- + REDUCTION OF PAIN SENSITIVITY**  
20% nitrous oxide has the same analgesic potency as 15mg morphine sc.
- + A MORE GENTLE SEDATION**  
No anesthesia – patient is conscious and cooperative.
- + REDUCED GAGGING REFLEX**  
Decreased gagging reflex and retching makes treatment faster and more comfortable.
- + EFFICIENCY & TIME SAVING**  
A patient with less anxiety is more cooperative and can be treated much more efficiently.
- + TITRATION**  
Optimal sedation by individually adjusted nitrous oxide mixture!
- + ALMOST NO RISKS**  
 $N_2O$  is a sedative with the lowest side effects!
- + NO ADVERSE EFFECTS**  
Gentle sedation. No metabolism (< 0.1 %).  
Recovery occurs within a few minutes.  
Fit to drive within 15 minutes.
- + MILLIONS OF PATIENTS**  
Millions of patients worldwide have been treated with nitrous oxide and have overcome their anxieties regarding dental treatment. Almost everyone would choose nitrous oxide again.
- + MARKETING & SERVICE**  
Nitrous Oxide Sedation can be a unique selling point with a positive marketing effect for your practice.

\* Reference: Müller, Ingrid: Angst essen Zähne auf.  
In: Focus-Gesundheit September/October 2015



  
**US Patented**  
No. 10,065,159 B2



Tested and approved in Germany!

# OXYGEN-NITROUS OXIDE MIXER

» The core of each Nitrous Oxide Sedation, the Mixer Unit.

Baldus presents two types of mixer units from its own production. The Baldus® Analog with its elegant design and its intuitive handling. The digital mixer unit Baldus® Touch combines the latest technology and innovative features and is a world debut on the dental nitrous-oxide market throughout the world.

What is the difference between the Baldus mixer units and other Oxygen - nitrous oxide mixer units?

**Created in Germany**  
**Designed in Germany**  
**Made in Germany**

## The world debut Baldus® Touch

### Timeless elegance combined with high technology – the Baldus® Touch.

Our everyday lives are unimaginable without the use of touchscreens. The Baldus® Touch is the next generation of digital oxygen - nitrous oxide mixer units. The new display menu functions and system settings give the Baldus® Touch the status of a modern and distinctive mixer unit on the dental nitrous-oxide market.

#### + US PATENTED

#### + HIGH-TECHNOLOGY

- **Touch-function** with a brilliant screen resolution.
- **Easy documentation** for the patient's treatment parameters and patient's ID.
- **Thermic measurement** works essentially independent from pressure and temperature, so therefore guaranties a shorter reaction time.
- **Software & Hardware Watchdog** guaranty safety.

#### + BALDUS® TOUCH FEATURES

- **Recovery button:** at the end of all sedation treatments the recovery button is pushed, which will produce a flow of 100 percent oxygen and a countdown will start.
- **Error memory:** in the documentation memory there is an error memory, which records the last 20 erros, which can be used for an analysis with the dentist and the Baldus service team.
- **Maximum total flow of 18 l/min:** The breathing bag will always be filled sufficiently, especially for patients with anxiety problems or patients with a bigger lung volume.

#### + HIGH SAFETY STANDARDS

- **Self-test** at power on, also in Standby mode.
- **Nitrous oxide lock** when lack of oxygen.
- **Acoustic and visual information signals.**
- **User-friendly single button control.**
- **Red O<sub>2</sub>-Flush-Button** for immediate, 100% oxygen supply.
- **Non-rebreathing valve** on the Bag-T.
- **Ambient air valve** on the Bag-T.
- **O<sub>2</sub> Demand Mask-Port.**
- **Standard coding of the gas types.**

## » The 21<sup>st</sup> Century Mixer Unit The Baldus® Touch



#### + HYGIENIC DESIGN

- **Fast and easy cleaning** due to flat fronts and rounded corners.
- **Intelligent placement** of the ambient air valve under the Bag-T protects from dust and other pollutants.
- **Filters** where the gas enters to keep any impurities away.
- **Microbiological testing:** requirements regarding the hygienic standards during production and packing are 100% fulfilled.





- + US Patented
- + Touchscreen
- + Automatic data documentation
- + Automatic lead out by recovery button
- + Maximum total flow of 18 l/min
- + O<sub>2</sub>-Flush-Button
- + Self-test-function & error memory
- + Software & Hardware Watchdog
- + Thermic measurement



### Mixer Unit „Baldus® Touch O<sub>2</sub>/N<sub>2</sub>O Mixer for Analgesia“

**Details:**

- Digital mixer unit controlled via HQ-touchscreen
- Max. mixing concentration 50%, 60% or 70% N<sub>2</sub>O (Standard delivery 70%)
- Mixing ratio is precisely settable over the % arrow key
- Automatic adjustment of the outcoming flow (user-friendly single-button control)
- Oxygen Flush for an immediate oxygen supply over a red-signal button
- Upgraded premium safety package with acoustical and visual information signals
- Self-test-function
- Automatic nitrous-oxide stop in case of lack of oxygen
- Non-rebreathing valve on the Bag-T
- Ambient air valve (when reservoir bag empty)
- Connection coding of the gas types
- O<sub>2</sub>-Demand Mask-Port
- USB-Port
- Oxygen nitrous oxide mixer unit offers the possibility of installation in the All-in-One-Cart 15005T, in a hollow wall or in the door of a working unit
- Input voltage: 110/ 230V AC, 50-60Hz
- Dimensions (Height x Width x Depth): 136 x 300 x 119 mm
- Weight: 2,5 kg
- US Patented No. 10,065,159 B2
- German Innovation Award Winner 2019

Delivery contains Mixer Unit and Bag-T; without mounting and scavenger system.  
**3 years warranty when concluding a maintenance agreement.**

Art.-No.	Article
16010-50	Baldus® Touch-50% N <sub>2</sub> O
16010-60	Baldus® Touch-60% N <sub>2</sub> O
16010-70	Baldus® Touch-70% N <sub>2</sub> O

**4.495,00 €**

## Baldus® Analog

### + US PATENTED

### + HIGH-QUALITY ROBUST MATERIAL

- Rugged in design and with over 92% of the mixer unit manufactured in **aluminium and stainless steel**

### + HYGIENE

- More **hygienic and easy to clean** due to flat surfaces with less prominent edges
- **The installation comes with a fresh air valve** which is under the Bag-T to prevent infiltration of dust and other pollutants

### + NEW FEATURES

- **Acoustic signal** in the event of a lack of oxygen (decrease in pressure under 3 bar)
- **Tamper-resistant:** Opening the mixer is only possible with a purpose-built-in tool
- Due to a **maximum flow** of 15 l/min the bag is always sufficiently filled
- **Centrally-placed an O<sub>2</sub> flush button** for a 100% oxygen supply and a quick reaction
- **Easily-readable nitrous oxide concentration:** Baldus® Analog displays the percentage of nitrous oxide concentration and not the oxygen concentration → no need to convert the gases concentration, therefore less misjudgement



### Mixer Unit „Baldus® Analog O<sub>2</sub>/N<sub>2</sub>O Mixer for Analgesia“

#### Details:

- Max. mixing concentration 50%, 60% or 70% N<sub>2</sub>O (Standard delivery: 70% N<sub>2</sub>O)
- CE-certified
- Mixing concentration precisely adjustable by %-N<sub>2</sub>O regulator
- Automatic adjustment of output flow in case of increased N<sub>2</sub>O-intensity (single-button operation)
- ISO colour coding
- Red marked O<sub>2</sub> flush button at the front of the device for immediate oxygen supply
- Premium safety package
- Automatic nitrous oxide stop and advisory tone in case of a lack of oxygen
- Non-rebreathing valve
- Ambient air valve (when the reservoir bag empty)
- Coding of the gas types for connection
- O<sub>2</sub> Demand Mask-Port
- Connection for emergency ventilation mask
- Width x Length x Height: 155 x 148 x 270 mm
- US Patented No. 10,065,159 B2

Delivery contains Mixer Unit and Bag-T; without mounting and scavenger system.  
**3 years warranty when concluding a maintenance agreement.**

Art.-No.	Article
16009-50	Baldus® Analog-50% N <sub>2</sub> O
16009-60	Baldus® Analog-60% N <sub>2</sub> O
16009-70	Baldus® Analog-70% N <sub>2</sub> O

**2.995,00 €**



*Work smarter.  
Not harder.  
Baldus® Lachgassedierung.*



**MADE  
IN  
GERMANY**

» OUR TIP

Do not forget oxygen and nitrous oxide to order.  
All information and prices see page 30.

## COMPREHENSIVE SET NITROUS OXIDE SEDATION FOR ADULTS

For a smooth start with Nitrous Oxide Sedation we offer a comprehensive set. We provide an overview about the important elements of the Nitrous Oxide Sedation that will help you start successfully in your practice.

### The package price includes

Item no.	Product	Baldus® Touch Set	Baldus® Analog Set
16010	Mixer Unit „Baldus® Touch O <sub>2</sub> /N <sub>2</sub> O Mixer for Analgesia“	✓	
16009	Mixer Unit „Baldus® Analog O <sub>2</sub> /N <sub>2</sub> O Mixer for Analgesia“		✓
16024-1	Baldus® Adult Double Mask Scavenger System, the package includes: autoclavable outer mask and inner mask, complete hose system and hose connectors	✓	✓
16002	Baldus® Vacuum Control Block with regulator valve	✓	✓
15015-3L	3-Liter Adult Reservoir Bag, latex free and autoclavable	✓	✓
15020	Ultraslim Cylinder Pressure Regulator O <sub>2</sub> including medical hose	✓	✓
15021	Ultraslim Cylinder Pressure Regulator N <sub>2</sub> O including medical hose	✓	✓
15004S	All-in-One-Cart Ultra Slim, encased portable cart for Baldus® Touch and Analog O <sub>2</sub> /N <sub>2</sub> O Mixer, Dimensions: 108 x 38,5 x 29 cm, 25 kg, Colour: matt white	✓	✓
	Marketing-Package: video with english subtitles for waiting room, information flyer for patients and medical devices book	✓	✓
<b>Package Price</b>		<b>7.395,- €</b>	<b>5.895,-</b>

# Comprehensive Set Baldus® Touch All-in-One Ultra-Slim



Package Price  
Baldus® Touch  
**7395,- €**  
Net price plus VAT/taxes

# Comprehensive Set Baldus® Analog All-in-One Ultra-Slim



Package Price  
Baldus® Analog  
**5.895,- €**  
Net price plus VAT/taxes



# BALDUS® SCAVENGER BREATHING SYSTEM

» The eco-friendly companion for the Nitrous Oxide Sedation.

The double nasal mask system from Baldus® is the „green“ innovation on the world nitrous oxide market. The development and production is 100% made in Germany.



## What makes the Baldus® double nasal mask scavenger system one of a kind?

### + SUSTAINABILITY

- Eco-friendly mask system.
- All parts – the tube system and the outer and inner mask with all the connectors – are autoclavable.
- Cost-saving – no consumption items.

### + HANDLING

- Unrestricted space for the dentist and his team.
- Easy handling.
- Markings and different sizes for an easy set up procedure.
- The patient's breathing can be constantly observed due to the translucent materials of the outer mask.
- The Baldus® Scavenger System is compatible with all oxygen-nitrous oxide mixer units.

### + SAFETY

- Special focus was put on the Scavenger System regarding leaking gases and air contamination during design development.
- The choice of userfriendly flexible materials ensures a perfect fit to all patient's face shapes.
- For greater safety the connectors have a bridge that avoids any mistakes when assembling the masks.

fully  
autoclavable



- ✓ Sustainable
- ✓ Eco-friendly
- ✓ High-quality

The Baldus® double nasal mask scavenger system – the eco-friendly system for the inhalation sedation with nitrous oxide.

**MADE  
IN  
GERMANY**



## + MATERIAL

- Use of translucent, soft and supple silicone.
- Latex free and autoclavable.
- Biocompatibility confirmed.

## + PATIENT COMFORT

- The userfriendly material prevents uncomfortable pressure marks on the skin of the patient.
- An oval membrane allows a minimal resistance during exhaling, which ensures a better and easy discharge exhaled gases.
- Thanks to the bigger volume of the inner masks, which allows optimal flow for the patient's inhaling and exhaling.


### Operating mode of the mask system:

How does the Baldus® Scavenger Breathing System work? Because of the exhalation in the inner mask the membrane opens and the exhaled gas enters the outer mask. The suction in the outer mask is separate from the inner mask, so the patient can breathe without resistance. This reduces the straw effect from other mask systems.

The double nasal mask system protects the dentist and his team successfully from the exhaled gas. The supple, autoclavable nasal masks are easy to change and available for adults and children, which are available in different colours and scents: Blueberry Dance, Bouncy Bubble, Happy Apple, Strawberry Delight and Chilla Vanilla.




## Baldus® Double Mask Scavenger System

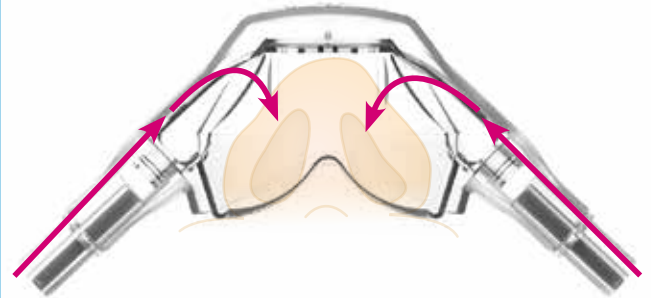


### Inhalation Phase

Delivery of the oxygen-nitrous-mixture.

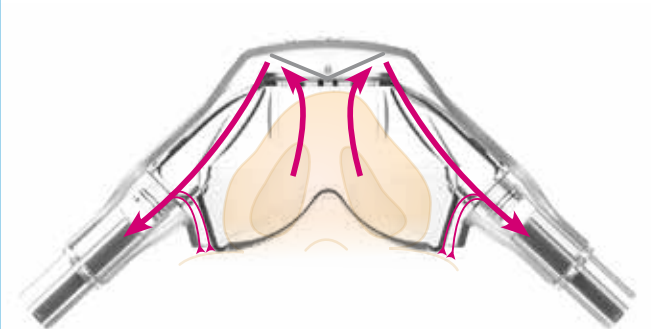

### Exhalation Phase

The exhaled gas flows via the oval membrane and the aerodynamic located channels optimal to the vacuum channel.

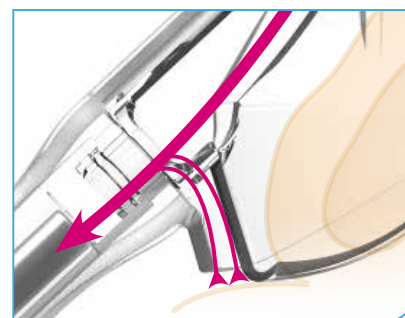

### Exhalation Phase

Minimized loss of N<sub>2</sub>O molecules because of the accurate fitting. The few molecules that leave the sealed Baldus® inner mask are sucked off from the unique deep-fitted Baldus® outer mask via the vacuum.

» BASIC-PACKET



**Baldus® Double Mask Scavenger System (CE-certified)**

**Details:**

- Autoclavable outer and inner mask, vacuum hose, double hose, medical gas supply hose and hose connectors

Art.-No.	Article
16024-1	Baldus® Double Mask Scavenger System, Adult
16030-1	Baldus® Double Mask Scavenger System, Pedo

550,00 €

» UPGRADE I



**Baldus® Double Mask Scavenger System + Vacuum Control Block (CE-certified)**

**Details:**

- Autoclavable outer and inner mask, double hose, medical gas supply hose and hose connectors
- Vacuum control block for controlling and adjusting the suction, including vacuum hose and adapter for saliva ejector

Art.-No.	Article
16024-2	Baldus® Double Mask Scavenger System + Vacuum Control Block , Adult
16030-2	Baldus® Double Mask Scavenger System + Vacuum Control Block , Pedo

699,00 €

**Material:** Nasal Masks, Hose Systems: medical silicone; Vacuum control block, Connectors: synthetic material  
**Compatible for devices of** Baldus®, Matrx, Porter, Accutron, Tecno Gaz



# The Nitrous Oxide Sedation THE HIGHLIGHT IN YOUR PRACTICE

The perfect education for patients with Sophie and Dr. Hansen.

» This is our free of charge service for your waiting room or website with a clear explanation of how nitrous oxide treatment works. Just download the explanatory video or contact us.



# Autoclavable Baldus® Double Nasal Masks

Re-usable!

**MADE  
IN  
GERMANY**



## Baldus® Double Mask

**Details:**  
 • Autoclavable outer and inner mask and hose connectors

Art.-No.	Article
16025	Baldus® Double Mask Adult
16031	Baldus® Double Mask Pedo

**279,00 €**



## Baldus® Outer Mask

**Details:**  
 • Autoclavable outer mask with hose connectors for Baldus® Double Mask Scavenger System, without inner mask

Art.-No.	Article
16026	Baldus® Outer Mask Adult
16032	Baldus® Outer Mask Pedo

**185,00 €**



## Baldus® Inner Mask

**Details:**  
 • Autoclavable inner mask for Baldus® Double Mask Scavenger System, with membrane

Art.-No.	Article
16027	Baldus® Inner Mask Adult
16033	Baldus® Inner Mask Pedo

**99,00 €**

**Material:** Nasal Masks: medical silicone  
**Compatible for:** Double mask scavenger systems of Baldus®, Porter and Accutron; further systems on request.

Mixer Units

Breathing System

Equipment

Emergency

Central Gas Supply

Maintenance

n<sub>2</sub>O-Academy



## BALDUS® SINGLE USE INNER MASKS

**Baldus® stands for quality, service and innovation.**

Scented masks excite not only kids, but also the adults. The inner masks for our sensitive noses can enjoy the pleasant natural fragrances during the nitrous oxide sedation – which has no effect on the Baldus quality. Exact opposite – Baldus has managed to establish new standards: simplified handling for you and your team, better absorption of the N<sub>2</sub>O gas molecules and better wearing comfort for the patient. Absolutly latex free.



## Baldus® Inner Masks



### Baldus® Inner Masks

**Scent:** Happy Apple  
**Content:** 25 pieces per aroma bag  
 in decorative box  
**Application:** SINGLE USE

## Happy Apple

Art.-No.	Article
16028-AA	Baldus® Inner Masks Adult
16028-PA	Baldus® Inner Masks Pedo

**100,00 €**




### Baldus® Inner Masks

**Scent:** Blueberry Dance  
**Content:** 25 pieces per aroma bag  
 in decorative box  
**Application:** SINGLE USE

## Blueberry Dance

Art.-No.	Article
16028-ABlue	Baldus® Inner Masks Adult
16028-PBlue	Baldus® Inner Masks Pedo

**100,00 €**

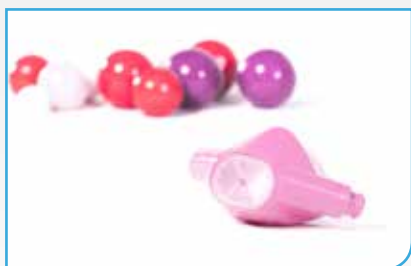

### Baldus® Inner Masks

**Scent:** Bouncy Bubble  
**Content:** 25 pieces per aroma bag  
 in decorative box  
**Application:** SINGLE USE

## BOUNCY BUBBLE

Art.-No.	Article
16028-ABub	Baldus® Inner Masks Adult
16028-PBub	Baldus® Inner Masks Pedo

**100,00 €**



### Baldus® Inner Masks

**Scent:** Strawberry Delight  
**Content:** 25 pieces per aroma bag  
 in decorative box  
**Application:** SINGLE USE

## Strawberry Delight

Art.-No.	Article
16028-AE	Baldus® Inner Masks Adult
16028-PE	Baldus® Inner Masks Pedo

**100,00 €**



### Baldus® Inner Masks

**Scent:** Chilla Vanilla  
**Content:** 25 pieces per aroma bag  
 in decorative box  
**Application:** SINGLE USE

## Chilla Vanilla

Art.-No.	Article
16028-AV	Baldus® Inner Masks Adult
16028-PV	Baldus® Inner Masks Pedo

**100,00 €**

**NEW**  
 unscented



### Baldus® Inner Masks

**Scent:** PURE GHOST - (unscented)  
**Content:** 25 pieces per aroma bag  
 in decorative box  
**Application:** SINGLE USE

## PURE GHOST

Art.-No.	Article
16028-AG	Baldus® Inner Masks Adult
16028-PG	Baldus® Inner Masks Pedo

**100,00 €**



## Baldus® Inner Masks

» In decorative acrylic glass

### Details:

- Application SINGLE USE
- 25 pieces per acrylic glass with comic label
- 25 buttons suitable to the scent
- High x diameter: 330 x 125 mm

## Baldus® LED Board

### Details:

- LED Board with the equipment for 3, 5 or 6 acrylic glasses
- Board colour: standard delivery in matt white RAL 9010
- Dimensions:  
3's Baldus® LED-Board: ca. 52,5 x 21,5 x 4 cm  
5's Baldus® LED-Board: ca. 65 x 21,5 x 4 cm  
6's Baldus® LED-Board: ca. 100 x 21,5 x 4 cm
- LED Light RGB with remote control
- Power supply 230V~ / 12V

Delivery without the acrylic glasses.

from **350,00 €**

## Baldus® Big Package

» Masks + LED Board

### Details:

- Baldus® LED Board
- 3, 5 or 6 filled acrylic glasses at the customer's choice (Art.-No. 16029)

from **699,00 €**

Art.-No.	Scent		Article	Price
16029-AA	Apple	<i>Happy Apple</i>	Baldus® Inner Masks Adult	129,00 €
16029-PA			Baldus® Inner Masks Pedo	129,00 €
16029-ABlue	Blueberry	<i>Blueberry Dance</i>	Baldus® Inner Masks Adult	129,00 €
16029-PBlue			Baldus® Inner Masks Pedo	129,00 €
16029-ABub	Bubblegum	<i>BOUNCY BUBBLE</i>	Baldus® Inner Masks Adult	129,00 €
16029-PBub			Baldus® Inner Masks Pedo	129,00 €
16029-AE	Strawberry	<i>Strawberry Delight</i>	Baldus® Inner Masks Adult	129,00 €
16029-PE			Baldus® Inner Masks Pedo	129,00 €
16029-AV	Vanilla	<i>Chilla Vanilla</i>	Baldus® Inner Masks Adult	129,00 €
16029-PV			Baldus® Inner Masks Pedo	129,00 €
16029-AG	unscented	<i>PURE GHOST</i>	Baldus® Inner Masks Adult	129,00 €
16029-PG			Baldus® Inner Masks Pedo	129,00 €



Size	Baldus® LED-Board		Baldus® Big Package	
	Art.-No.	Price	Art.-No.	Price
3's	16045	350,00 €	16045X	699,00 € 737,00 €
5's	16034	450,00 €	16034X	949,00 € 1.095,00 €
6's	16046	550,00 €	16046X	999,00 € 1.273,00 €

# BALDUS® MOUNTING SYSTEMS

MADE  
IN  
GERMANY

For mounting the Mixer Units we offer several options. If a cabinet installation or wall mount are not desired, we have the ideal mobile system for the oxygen-nitrous oxide mixer. The Baldus® All-in-One-Cart is perfect for every practice situation.



Available  
in all  
RAL colours



## All-in-One-Cart Ultra Slim

» Encased All-in-One-Cart for Mixer Baldus® Analog

### Details:

- Encased portable Cart with double solid rubber tyres for top driving comfort
- Mounting of the gas cylinders in the cart (gas cylinders are completely covered)
- Fulfilment of the safety requirements according to council directive 93/42/EEC
- Mounting of the Mixer Unit at the front of the cart
- Colour: matt white RAL 9010, further colours available upon request
- Height x Width x Depth: 108 x 38,5 x 29 cm
- Wight: 25 kg (without sedation system)
- Made in Germany

Art.-No.	Article
15004AS	All-in-One-Cart Baldus® Analog

Price without sedation system. **1.285,00 €**



## All-in-One-Cart Ultra Slim

» Encased All-in-One-Cart for Mixer Baldus® Touch

### Details:

- Encased portable Cart with double solid rubber tyres for top driving comfort
- Mounting of the gas cylinders in the cart (gas cylinders are completely covered)
- Fulfilment of the safety requirements according to council directive 93/42/EEC
- The Mixer Unit is integrated in the cart
- Colour: matt white RAL 9010, further colours available upon request
- Height x Width x Depth: 108 x 38,5 x 29 cm
- Wight: 25 kg (without sedation system)
- Made in Germany

Art.-No.	Article
15004TS	All-in-One-Cart Baldus® Touch

Price without sedation system. **1.285,00 €**

## » UPGRADES FOR ALL-IN-ONE-CART

Art.-No.	Article	Price
15010G	modern glass plate (float glass)	<b>195,00 €</b>
15010GRAL	modern glass plate, in all RAL colour standards	<b>240,00 €</b>
15010W	high-gloss finish, pure white	<b>600,00 €</b>
15010S	high-gloss finish, in all RAL colour standards	<b>680,00 €</b>



### All-in-One-Cart PIN INDEX

» Encased All-in-One-Cart for Mixer Baldus® Touch

#### Details:

- Space for **4 gas cylinders** (2 Oxygen & 2 Nitrous Oxide) **E-Size**
- Encased portable Cart with double solid rubber tyres for top driving comfort
- Mounting of the gas cylinders in the cart (gas cylinders are completely covered) CE-certified
- Fulfilment of the safety requirements according to council directive 93/42/EEC

- The Mixer Unit is integrated in the cart
- Standard Colour: matt white RAL 9010, further colours available upon request
- Made in Germany
- Upgrades page 23

Art.-No.	Article
15004T4	All-in-One-Cart

Price without sedation system. **1.285,00 €**



### Mobile stand compact 2-in-1

#### Details:

- Mobile stand is well suitable for the comfortable transportation of the analog mixer unit as well as the oxygen and nitrous oxide gas cylinders

- Free mobility of the mobile stand
- Fixable with a break
- Dimensions (H x W): 1400 x 630 cm
- Load: 60 kg

Art.-No.	Article
16041	Mobile stand compact 2-in-1

**779,00 €**



### Baldus® Mobile Cart

#### Details:

- Mobile mount
- Adjusts 360 mm - 560 mm
- 5-leg base
- CE-certified

- Compatible for Baldus® Analog, Matrix Quantiflex and Porter Flowmeter

Art.-No.	Article
16021	Baldus® Mobile Cart

**475,00 €**





## Baldus® Telescoping Swing Arm Bracket

### Details:

- Telescoping swing wall arm (180°)
- Plate measures: 200 x 75 mm
- Telescoping arm:  
adjusts 280 mm - 445 mm

- Including mounting peg
- Compatible for Baldus® Analog, Matrx Quantiflex and Porter Flowmeter

Art.-No.	Article
16020	Baldus® Telescoping Swing Arm Bracket

**378,00 €**



## Non-Telescoping Wall Bracket or Cart Mount

### Details:

- Mounting system for Mixer Baldus® Analog – wall or All-in-One-Cart
- Including fixing screws

Art.-No.	Article
16037	Wall Bracket or Cart Mount

**199,80 €**



## Bottle-Trolley

### Details:

- Bottle-Trolley for a comfortable transport of nitrous oxide and oxygen gas bottles
- Free movement of the Bottle-Trolley
- Height x Width x Depth: 97 x 60 x 50 cm
- Loading capacity up to 60 kg

Art.-No.	Article
15016	Bottle-Trolley

**299,00 €**



## Cylinder Wall Bracket

### Details:

- Wall Bracket with bottle clamp
- Holds for one or two gas bottles

Art.-No.	Article	Price
H547-1000	Wall Bracket for 1-Gas cylinder	<b>38,93 €</b>
H547-1100	Wall Bracket for 2-Gas cylinder	<b>49,70 €</b>

from **38,93 €**



## SEDATION EQUIPMENT

Our wish is that your patients feel safe in your good hands, which will always bring the patients back for further treatments. We offer you selective, high quality and certified equipment to complete your Nitrous Oxide Sedation System.



### Baldus® Vacuum Control Block

**Details:**

- CE-marked
- Vacuum control block with regulator valve for controlling and adjusting the suction of the Nitrous Oxide Sedation System (white/transparent)
- Installation between vacuum hoses
- Including autoclavable vacuum hose

Art.-No.	Article
16002	Baldus® Vacuum Control Block

**149,00 €**



### Reservoir Bag

**Details:**

- CE-marked
- Autoclavable and latex free
- Steam sterilization: 134°C, 5 min.

Art.-No.	Article
15015-3L	3-Liter Reservoir Bag (Adult)
15015-2,3L	2,3-Liter Reservoir Bag (Pedo)
15015-1L	1-Liter Reservoir Bag

**66,00 €**



### Nitrous Evacuator

**Details:**

- Active vacuum pump
- Nitrous oxide discharged directly into the open air
- Suitable for all scavenger systems
- Alternative of suction device on dental chair

- Height x Width x Depth: 227 x 302 x 154 mm
- Weight: 5,5 kg
- Voltage: 220-240V
- Frequency: 50/60Hz

Art.-No.	Article
15030	Nitrous Evacuator

**910,00 €**



## Ultraslim Pressure Reducer O<sub>2</sub>

» Hand connection for O<sub>2</sub> pressure cylinders of Nitrous Oxide Sedation System

### Details:

- Adjusted to 4,5 bar
- Including medical fabric hose
- DIN Cylinder Germany (NIST and DISS connection)
- Including adapter for local gas cylinder connection

Art.-No.	Article
15020	Ultraslim Pressure Reducer O <sub>2</sub>

**425,00 €**



## Ultraslim Pressure Reducer N<sub>2</sub>O

» Hand connection for N<sub>2</sub>O pressure cylinders of Nitrous Oxide Sedation System

### Details:

- Adjusted to 4,5 bar
- Including medical fabric hose
- DIN Cylinder Germany (NIST and DISS connection)
- Including adapter for local gas cylinder connection

Art.-No.	Article
15021	Ultraslim Pressure Reducer N <sub>2</sub> O

**425,00 €**



## Regulators with Connection for O<sub>2</sub> and N<sub>2</sub>O

» Regulators International Fittings:  
PIN INDEX, BULLNOSE, BENELUX BULLNOSE, AFNOR, UNI 4406

### Details Oxygen PIN-INDEX Pressure Reducer:

- Inlet, 1 x DISS Outlet
- Non Therapy

Art.-No.	Article
Divers	Regulators with Connection

**425,00 €**

### Details Nitrous Oxide PIN-INDEX Pressure Reducer:

- Inlet, 1 x DISS Outlet
- Non Therapy



## Medical Flow Meter

» Oxygen Medical Pressure Regulator

### Details:

- Content pressure gauge
- Flow rates: 0 – 1 – 2 – 3 – 4 – 5 – 6 – 7 – 8 – 9 – 12 – 15 L/min.
- 9/16" outlet
- Connection possibility for Bag-Valve-Mask Resuscitator, nasal cannula, etc.

- Including adapter for local gas cylinder connection

Art.-No.	Article
HDM01-15-KRPS	Medical Flow Meter

**460,00 €**

Alternative  
for Ultraslim  
Pressure  
Reducer O<sub>2</sub>



## Medical Flow Meter with DIN coupler

» Oxygen Medical Pressure Regulator

### Details:

- Content pressure gauge with separate outlet – DIN coupling 13260 (5,0 bar to 160 L/min.)
- 9/16" outlet
- Flow rates: 0 – 1 – 2 – 3 – 4 – 5 – 6 – 7 – 8 – 9 – 12 – 15 L/min.
- Connection possibility for Bag-Valve-Mask Resuscitator, nasal cannula, etc.

- Including adapter for local gas cylinder connection

Art.-No.	Article
AERO-way-Fast-15PS	Medical Flow Meter with DIN coupler

**470,00 €**



Image: Art.-No. 15006

## O<sub>2</sub> Pressure Hose

» Medical hose for the connection of O<sub>2</sub> pressure reducer or gas outlets to the mixer

### Details:

- Colour: ISO
- Length: max. 2,50 meters, other lengths available upon request
- Different connectors (see Article-Table)

Art.-No.	Article
15006	Medical hose – NIST and DISS connection (male)
16038-1	Medical hose – DIN gas probe and DISS connection
H6600107	Medical hose – DIN gas probe and NIST connection

89,00 €



Image: Art.-No. 16038-2

## N<sub>2</sub>O Pressure Hose

» Medical hose for the connection of N<sub>2</sub>O pressure reducer or gas outlets to the mixer

### Details:

- Colour: ISO
- Length: max. 2,50 meters, other lengths available upon request
- Different connectors (see Article-Table)

Art.-No.	Article
15007	Medical hose – NIST and DISS connection (male)
16038-2	Medical hose – DIN gas probe and DISS connection
H6600137	Medical hose – DIN gas probe and NIST connection

89,00 €



## O-Ring for pressure reducer

### Details:

- O-Ring for oxygen O<sub>2</sub> or nitrous oxide N<sub>2</sub>O standard pressure reducer
- Set of 10 O-Rings

Art.-No.	Article
1125714	O-Ring for oxygen O <sub>2</sub> standard pressure reducer
104746-4	O-Ring for nitrous oxide N <sub>2</sub> O standard pressure reducer

price per set: 38,00 €



## Nasal cannula flow set

### Details:

- Adapting connector and locking adapter for the Bag-T of the mixer unit
- 10 nasal cannulas

Art.-No.	Article
16040	Nasal cannula flow set

96,00 €



## Multiparameter Patient Monitoring Deluxe70

### Details:

- Oxygen saturation display, pulse rate, blood pressure, ECG, respiration and temperature
- TFT-Display in colour, 17,7 cm
- Suitable for hospitals, practices and private home
- For adults and children
- Free memory for 800 blood pressure readings and 24h ECG-curves
- Rechargeable battery
- Power connection possible
- 100-240VAC, 50/60Hz
- Visual and acoustic alarm
- Safety standard: IEC 60601-1
- Made in China

Art.-No.	Article
Deluxe70	Multiparameter Patient Monitoring Deluxe70

1.190,00 €



### Prince 100C1 Finger Clip Pulse Oximeter

**Details:**

- CE-certified: 0123
- Automatic ON/ OFF (finger sensor)
- Pulse rate and oxygen saturation display
- Plethysmogram display in real time
- Acoustical observable pulse
- Full-colour LED display
- Automatic screen rotation in 4 directions, easy to read
- Changeable to perfusion index display: display of finger circulation (low perfusion can lead to incorrect SPO2 signal)
- Visual and acoustic alarm
- Alarm limit individually adjusted
- Anti-motion technology
- Small and lightweight
- Low battery consumption (lasts up to 30 hours)
- Including 2 AAA batteries and device case
- Made in China

Art.-No.	Article
100C1	Fingerclip Adult
100D2	Fingerclip Child

**159,00 €**

Image similar



### SpO<sub>2</sub> Monitor (Pulsoximeter) with Finger Sensor OxyTrue A

**Details:**

- Oxygen saturation display (SpO<sub>2</sub>), pulse rate, plethysmogram, bar graph, short-term and long-term trend
- Alarm limit individually adjusted
- Intelligent linking of visual and acoustical information and signals
- Automatic suppression of surrounding light
- Signal quality display
- Fuzzy logic for extreme fast measurement display
- Intuitive operation with explicit error messages
- Data memory with integrated data management
- SpO<sub>2</sub> Monitor ideal for rescue service, device and sensor are extremely robust
- High speed data transmission via standard USB port
- Selection of 18 languages
- Use of device: continuous operation for more than 2 days / power saving mode up to 4 days
- Weight: approx. 160g (without batteries and sensor)
- Including SC 6500 Soft-Cap Sensor, PC software, USB data cable, silicone protector, batteries
- Made in Germany

Art.-No.	Article
15018	SpO <sub>2</sub> Monitor

**459,00 €**

**» UPGRADE FOR THIS ARTICLE**

Art.-No.	Article	Price
5825506	SoftCap Sensor, Pedo (10 to 20 kg)	<b>220,00 €</b>

# ADVERTISING PACKAGES

## The perfect basis for your own practice advertising!

### » BALDUS S



#### Patient Flyer

» Benefits of nitrous oxide sedation at a glance.


##### Details:

- Format: DIN A5
- Package: 100 flyer
- Individual adaptation of your practice CI upon request

Art.-No.	Article
16008	Baldus S Patient Flyer

**15,00 €**

### » BALDUS M

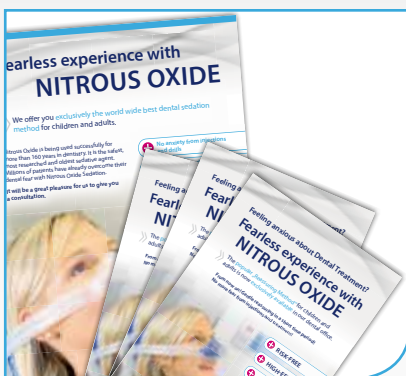

#### 3 Advertising Posters + 100 Patient Flyer

» Effective advertising for your waiting room.

##### Details:

##### Advertising Posters

- Format: DIN A2
- Package: 3 poster

##### Patient Flyer

- Format: DIN A5
- Package: 100 flyer
- Individual adaptation of your practice CI upon request

Art.-No.	Article
16007	Baldus M Advertising Package

**26,50 €**

### » BALDUS BUTTONS



#### 50 „Medals of Honor“

» The nitrous-oxide heroes as a individual medals for every kid.

##### Details:

- Format: 25 mm needle button
- Package: 3 poster

Art.-No.	Article
16040-A	Baldus Buttons „Happy Apple“
16040-Blue	Baldus Buttons „Blueberry Dance“

Art.-No.	Article
16040-Bub	Baldus Buttons „Bouncy Bubble“
16040-E	Baldus Buttons „Strawberry Delight“
16040-V	Baldus Buttons „Chilla Vanilla“
16040-G	Baldus Buttons „Pure Ghost“

**19,95 €**

### » BALDUS KIDS PREMIUM



#### 25 Colouring Books + 25 „Medals of Honor“

#### + 5 soft toy „Strawberry Delight“

» Create enthusiasm to every kid in 3 steps!

##### Details:

##### Before treatment:

- 25 Baldus® Colouring Books with coloured pencils (set 4 pcs).  
In the waiting room the children get to know their nitrous-oxide heroes.

##### After treatment:



- 25 Buttons. As a reward the children receive a button with their own nitrous-oxide hero.
- 5 soft toys „Strawberry Delight“ for very brave kids.

Art.-No.	Article
16039	Baldus Kids Premium

When placing an order, please specify the nitrous-oxide hero for the buttons.

**39,50 €**

# SCAVENGING CIRCUIT COMPONENTS

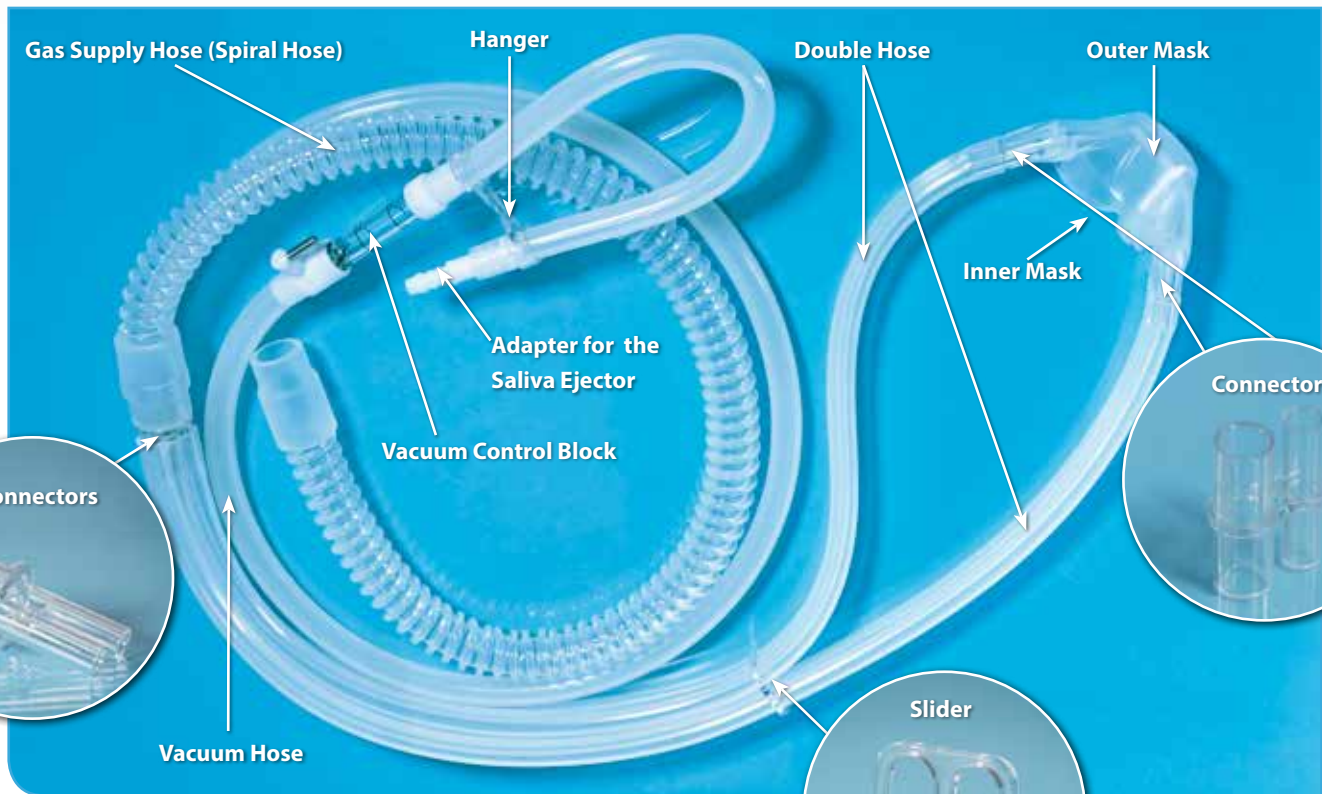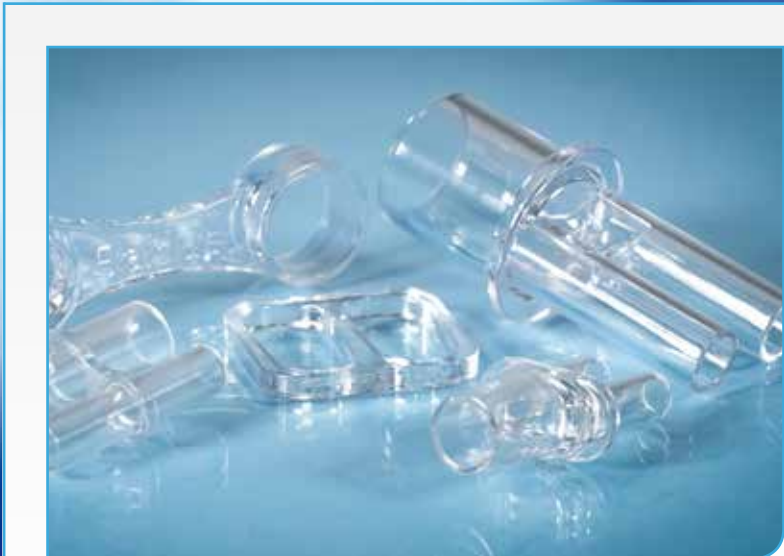

Mixer Units

Breathing System

Equipment

Emergency

Central Gas Supply

Maintenance

n<sub>2</sub>O-Academy



## Oval Membrane for Baldus® Masks

### Details:

- Content: 10 pieces
- Material: medical silicone, autoclavable
- Adapted for autoclavable and disposable inner masks

Art.-No.	Article
16036	Oval Membrane

price per set: **24,50 €**



## Baldus® Adapter for Saliva Ejector

### Details:

- White adapter for the saliva ejector of the small suction cup to the vacuum hose

Art.-No.	Article
16014	Baldus® Adapter for small Saliva Ejector

**7,85 €**



## Baldus® Y-Adapter for Gas Supply Hose (Spiral Hose)

### Details:

- Adapter for gas supply hose 22 mm to double hose
- Autoclavable

Art.-No.	Article
16016	Baldus® Y-Adapter for Gas Supply Hose

**26,80 €**



## Baldus® Y-Connector for Vacuum Hose

### Details:

- Adapter for vacuum hose 13,3 mm to double hose
- Autoclavable

Art.-No.	Article
16017	Baldus® Y-Connector for Vacuum Hose

**24,60 €**



## Baldus® Mask-Connectors

### Details:

- Connectors for double nasal mask to double hose
- Autoclavable
- Content: 2 pieces

Art.-No.	Article
16018	Baldus® Mask-Connectors

**27,35 €**





## Baldus® Vacuum Hose Hanger

### Details:

- Hanger for vacuum hose to read the vacuum control block vertically
- Autoclavable

Art.-No.	Article
16015	Baldus® Vacuum Hose Hanger

**29,80 €**



## Silicone-Vacuum Hose

### Details:

- Autoclavable
- Length: 2,5 Meters
- Compatible for Baldus® Scavenger Breathing System and Baldus® Vacuum Control Block

Art.-No.	Article
16003	Silicone-Vacuum Hose

**69,00 €**



## Baldus® Slider

### Details:

- Function: for fixing the Baldus® Scavenger Breathing System to the patient's head
- Autoclavable

Art.-No.	Article
16019	Baldus® Slider

**19,50 €**



## Silicone-Double Hose

### Details:

- Autoclavable
- Content: 2 Double Hoses
- Compatible for Baldus®, Porter and Accutron scavenger systems
- Length: 70 cm

Art.-No.	Article
16004	Silicone-Double Hose

**196,00 €**



## 22 mm Gas Supply Hose (Spiral Hose)

### Details:

- Autoclavable
- Length: 90 cm, 180 cm
- Compatible for Baldus® Scavenger Breathing System

Art.-No.	Article	Price
16013-90	90 cm Gas Supply Hose	<b>76,00 €</b>
16013-180	180 cm Gas Supply Hose	<b>89,00 €</b>

**from 76,00 €**




## Emergency Equipment

### FIRST AID

#### Saving Lives!

An emergency could occur unexpectedly – at any time and in any place. Everybody of all age groups can be affected; it does not matter if they are athletes, managers or students.

We offer high quality equipment for an optimal and saving emergency action.

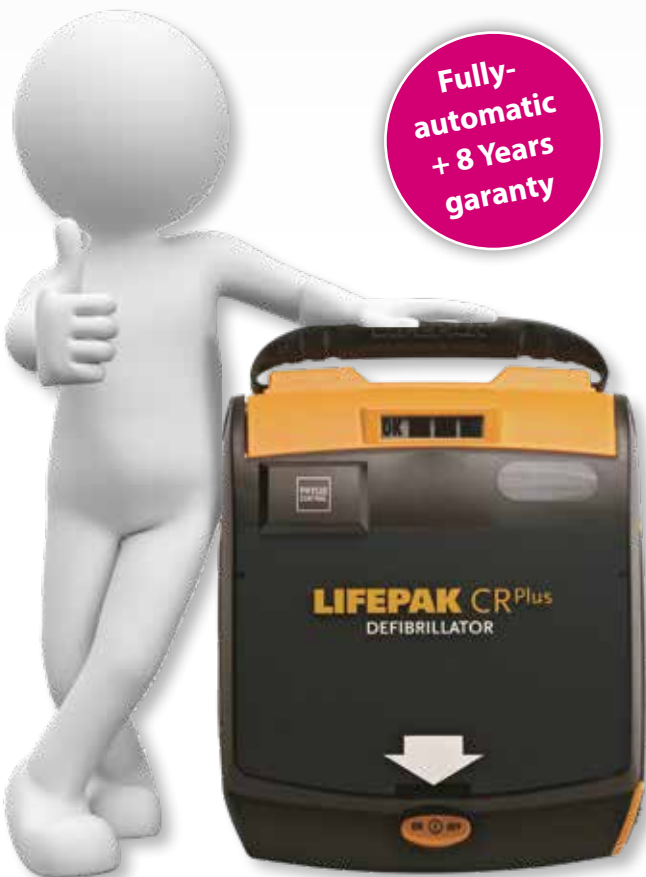


#### LIFEPAK CR® Plus AED fully-automatic

##### Details:

- Compact and very easy to use
- Automatic external defibrillator, English
- Simply open the cover, follow the verbal instructions and push the button to deliver a shock if necessary
- Loud voice prompts provide guidance on applying electrodes and initiating a shock
- Includes: one battery, 2 pairs electrodes, AMBU® First Responder Kit, carry case
- Weekly self-tests
- AED assembled and ready to use

Art.-No.	Article
201000	LIFEPAK CR® Plus AED

**1.399,00 €**

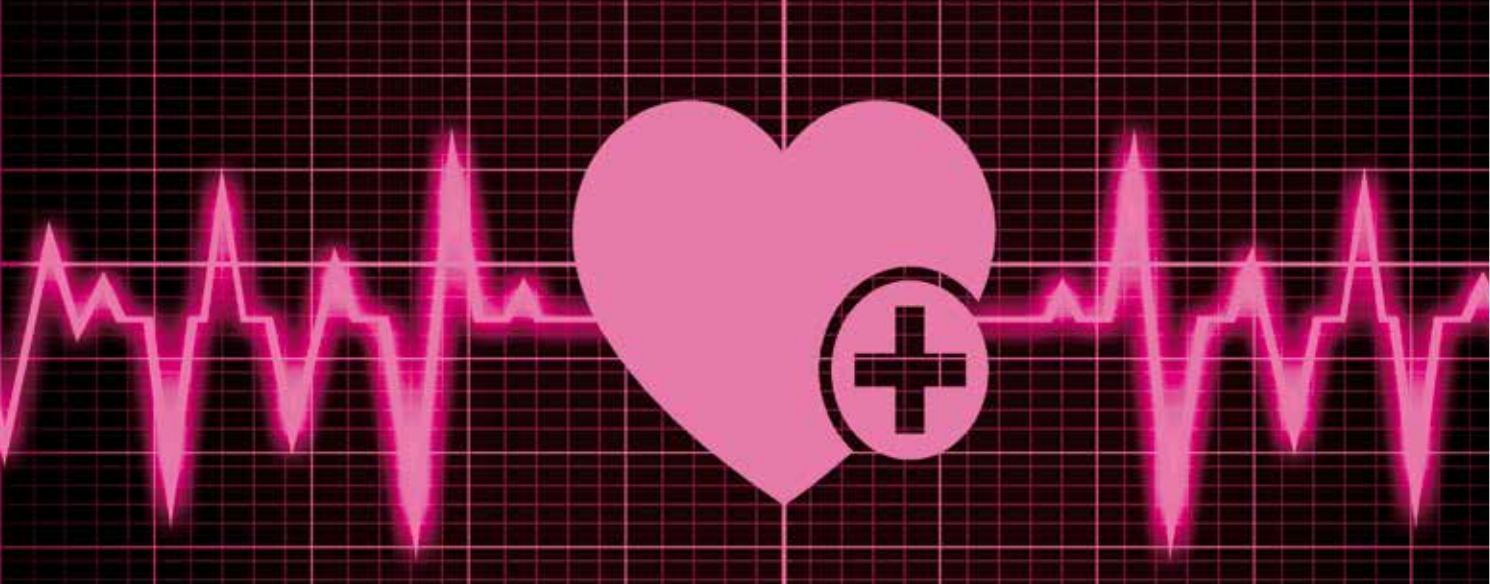

For your  
Nitrous Oxide  
Sedation



### Baldus® Demand Valve

**Details:**

- Demand Valve with universal mask, 3 meters hose and DISS connection for Baldus® Nitrous Oxide Sedation
- Ventilation active and passive

Art.-No.	Article
15012	Baldus® Demand Valve

**350,00 €**



### Disposable Bag-Valve-Mask Resuscitator

**Details:**

- Single Use
- Available in two models (Adult and Child)
- Including 1 mask

Art.-No.	Article
HBB06-E-45R	Adult, Disposable Mask Resuscitator
HBB06-K-23R	Child, Disposable Mask Resuscitator

**24,99 €**



### Disposable Ventilation Mask

**Details:**

- Available in sizes 1-5
- 4 masks per size

Art.-No.	Article
HBM06-UV	Disposable Ventilation Mask

Price per set: **12,00 €**



### Oxygen Tubing with Nasal Cannula

**Details:**

- Single Use
- Length 2,10 meters
- Adult size
- Latex free
- Supplied in packs of 10

Art.-No.	Article
14090510	Oxygen Tubing with Nasal Cannula

Price per Set: **25,00 €**



## Central Gas Supply

Our quality management and our ISO 13485 certification make it possible to carry out the planning, montage and service of medical gas supply systems.

**A central gas supply system offers several major advantages:**

- Cost saving – gas cylinders at a central location result in less handling
- Higher safety level because of high pressure gas storage and installation located in specified and safe place
- Uninterrupted gas supply with automatic changeover
- Warning signal at low gas pressure
- Visual integration in the modern design of your practice



Please contact us at any time. We will be pleased to support you with our expert's advice about the central gas supply and the possibilities to install the system in your practice.



**MADE  
IN  
GERMANY**



Mixer Units

Breathing System

Equipment

Emergency

Central Gas Supply

Maintenance

n<sub>2</sub>O-Academy



## Gas Manifold – Changeover System

» Fully automatic changeover manifold oxygen or nitrous oxide, AGS 35

### Details

- Changeover unit for regulation, supervision and control of medical and technical gas supply units
  - 2 x 1 cylinder with pneumatically-controlled change-over
  - 2 high capacity gas pressure regulators
  - Outlet nominal pressure 3-10 bar
  - 2 connection bends made of copper (copper pipe 8 x 1,5 mm)
  - Procedural message
  - Bottle bracket for 2 gas cylinders of 50 l volume
  - Signs-Set
- Further information available upon request.**




## Gas Manifold – Single Cylinder System

» Pressure reducer for one cylinder oxygen or nitrous oxide

### Details

- Outlet nominal pressure 3-10 bar
  - 1 connection bend made of copper (copper pipe 8 x 1,5 mm)
  - Procedural message
  - Bottle bracket for 1 gas cylinder of 50 l volume
  - Signs-Set
- Further information available upon request.**

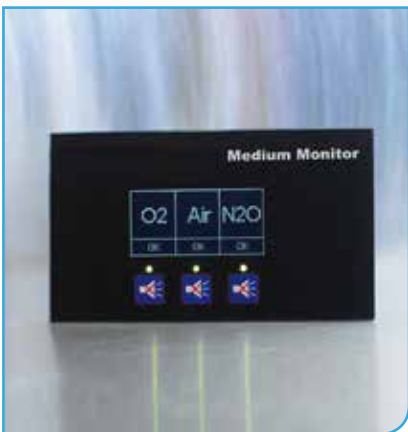


## Controllpanel for Manifold

» Medium Monitor equipped with an innovative glass display

### Details

- Indication via illuminated LCD display
  - Information on current status (OK and warning of high and low pressure)
  - Operation by means of touch keys
  - Visual and acoustic signals
  - Communication of several emergency signals via bus system
  - 3 potential-free contacts
- Further information available upon request.**



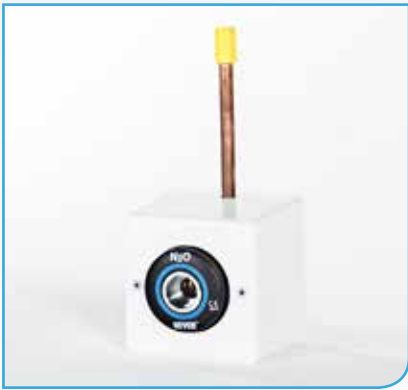

## Medical Gas Outlet O<sub>2</sub> / N<sub>2</sub>O

### Details

- Simple installation – for concealed installation, for surface (on plaster) installation
- Fast connection and disconnection
- Designed for medical environment, small size and easy to clean

- Complies with colour coding and description by standard
- DIN-Standard

Further information available upon request.



## Flowmeter DIN Direct

### Details

- Flow rates adjustable in 9 levels from 0 to 15 L/min.
- Supply pressure max. 5 bar
- Body made of anodized aluminium
- DIN Connector
- DISS O<sub>2</sub>-Outlet

- For direct connection to medical gas outlet
- Exit for nasal cannula or resuscitation bag

Art.-No.	Article
H6600400	Pressure Reducer Direct

**171,00 €**



24 CHROM-VAN





## Maintenance

### » Baldus Medizintechnik GmbH offers a maintenance and inspection service.

In order to protect both patients and staff, you have maintenance and safety checks carried out by Baldus Medizintechnik GmbH or an authorised service partner once a year.

The technical personnel of Baldus Medizintechnik travels throughout Europe for the highest priority in your practice is the medical care of the patients.

#### Maintenance and safety inspections of oxygen-nitrous oxide mixer

In the need of maintenance or repair, send in your oxygen-nitrous oxide mixer for inspection to Baldus Medizintechnik GmbH. This service is valid worldwide. For Europe our qualified personnel is also available for a visit in your practise.

Maintenance interval and safety inspection of oxygen-nitrous oxide mixer: annual service

#### Maintenance of central gas supply

We offer maintenances for the following central systems:

- Oxygen
- Nitrous Oxide
- Compressors, compressed air

A regular maintenance of the central gas supply is indispensable.

A continuous supply of medical gases is necessary for all patients. In case of a gas system breakdown it can have life-threatening consequences for patients and also be very costly.

#### Additional maintenances

The technical team of Baldus Medizintechnik carries out maintenance also for the following sectors:

- **Quality control of compressed air required by European Pharmacopoeia**
- **Terminal Units (gas outlets) Air, N<sub>2</sub>O, O<sub>2</sub>**
- **Anaesthetic scavenging systems**
- **Ceiling-mounted supply units**
- **Compressed air preparation for breathing air**

We control and maintain devices of the following manufacturers:

- Baldus®
- Porter
- Matrix
- Further producers on request

**Further information on request.**





## n<sub>2</sub>o-Academy Workshops

### TRAINING FOR NITROUS OXIDE SEDATION

» The professional use of nitrous oxide sedation. We will make you fit in the shortest time!

**Certification according to the Council of European Dentists resolution 2012.**

The n<sub>2</sub>o-Academy workshops include the presentation and explanation of all important aspects of nitrous oxide sedation. To be allowed to sedate with nitrous oxide, a dentist has to attend 10-12 hours workshop (CED resolution 2012).

#### Workshop Content:

- + Pain and anxiety in dental medicine
- + History of nitrous oxide
- + Characteristics and effects of nitrous oxide
- + Range of applications
- + Medical indications and contraindications
- + Legal aspects
- + Ambient air contamination and compliance
- + Equipment technology
- + Dental office marketing
- + Economic efficiency calculation; Billing
- + Role play & practical application
- + Emergency- and quality- management
- + Examination and Certificate handover

Our seminar instructors and lecturers are highly qualified experts with years of experiences with nitrous oxide sedation.

#### Workshop information:

For our workshops, we see a big importance in extensive theoretical information, but also the practical part!

Each member of your team has the opportunity to attend our workshop. Your assistant will be trained and qualified to apply the sedation method professionally and with high standards.

#### Private training worldwide:

Would you like to certify more employees and colleagues in the nitrous oxide sedation at your practice? No problem. We can offer our workshop as well as Inhouse workshop for up to 20 participants.

**Further information available on request!**

#### » INFORMATION

According to the Council of European Dentists every dentist has to visit a 2-day workshop (at least 10 hours) before using nitrous oxide for sedation.

## Trust the specialists in **MEDICAL GASES AND INHALATION SEDATION**

» Our service has been evaluated by doctors who use nitrous oxide as being a grade 1 (1.4), the highest possible grade

### » INFORMATION

All prices are net price, excluding shipment and TAX. Please pay attention to our regular sales conditions.

Baldus Medizintechnik GmbH will not be liable for patient treatment and improper handling. The treating physician is responsible for patient monitoring and needs to provide a certificate of competence. Please read the instructions for use carefully and in full before first operating the device in order to protect you and your patients from any incorrect operation.

Please note, that the images may deviate slightly from the products supplied.  
All information without engagement.

